

SHELL WIDTH & LENGTH DATA SHEET

ortesi

3D printed braces & orthotics

WIDTH



NARROW



STANDARD



WIDE

- NARROW MEASURED FROM THE BISECTION OF THE 5TH MET HEAD TO THE LATERAL ASPECT OF THE 1ST MET HEAD
- STANDARD MEASURED FROM THE LATERAL ASPECT OF THE 5TH MET HEAD TO THE BISECTION OF THE 1ST MET HEAD
- WIDE MEASURED FROM THE LATERAL ASPECT OF THE FOOT TO THE MEDIAL ASPECT OF THE 1ST MET HEAD


LENGTH



SHORT



STANDARD

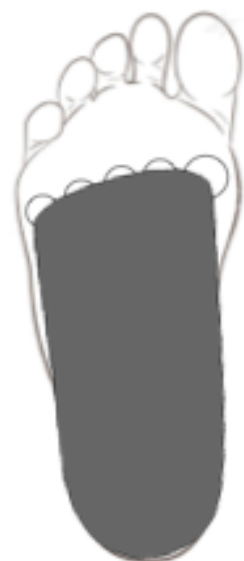

LONG

- SHORT MEASURED PROXIMAL TO THE 1ST AND 5TH MET HEADS
- STANDARD MEASURED PROXIMAL TO THE BISECTION OF THE 1ST AND 5TH MET HEADS
- LONG MEASURED DISTAL TO THE 1ST AND 5TH MET HEADS